

ATHLETIS™ Ultra-High-Pressure Balloon

ULTRA-HIGH PRESSURE. ULTRA-LOW PROFILE.



Optimized Deliverability to Reach and Treat Challenging Lesions



ULTRA-HIGH PRESSURE

Athletis' patented Vectran™ braided balloon design enables optimized deliverability and ultra-high pressures up to 40 ATM. The unique, non-compliant balloon design enables Athletis to open to the stated diameter even in highly resistant lesions.

Vectran™ is a high-performance material that was used by NASA in the Mars Rover parachute.

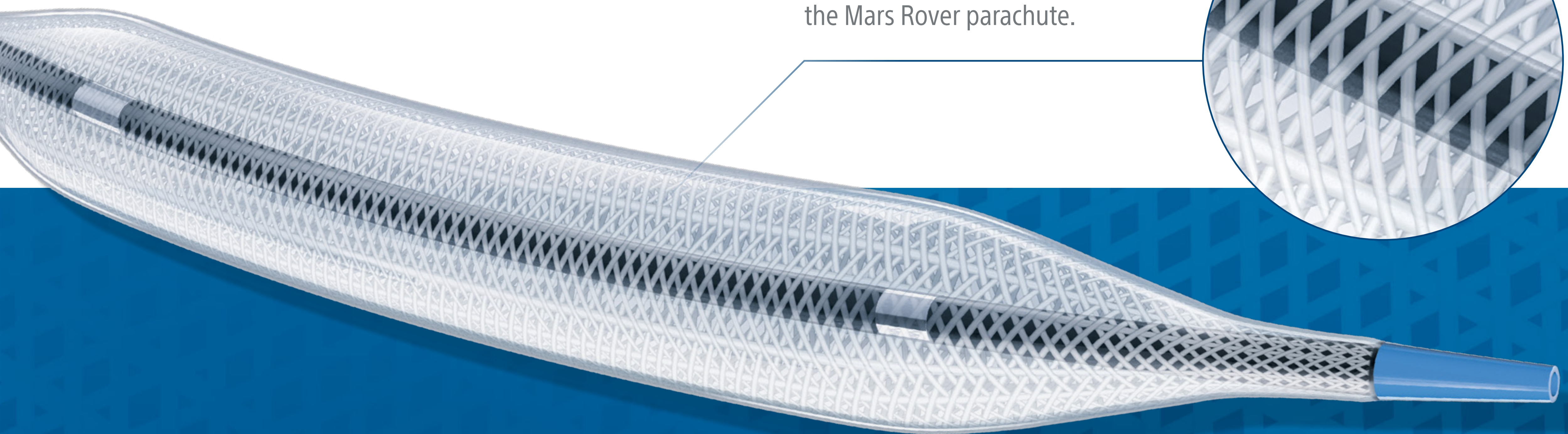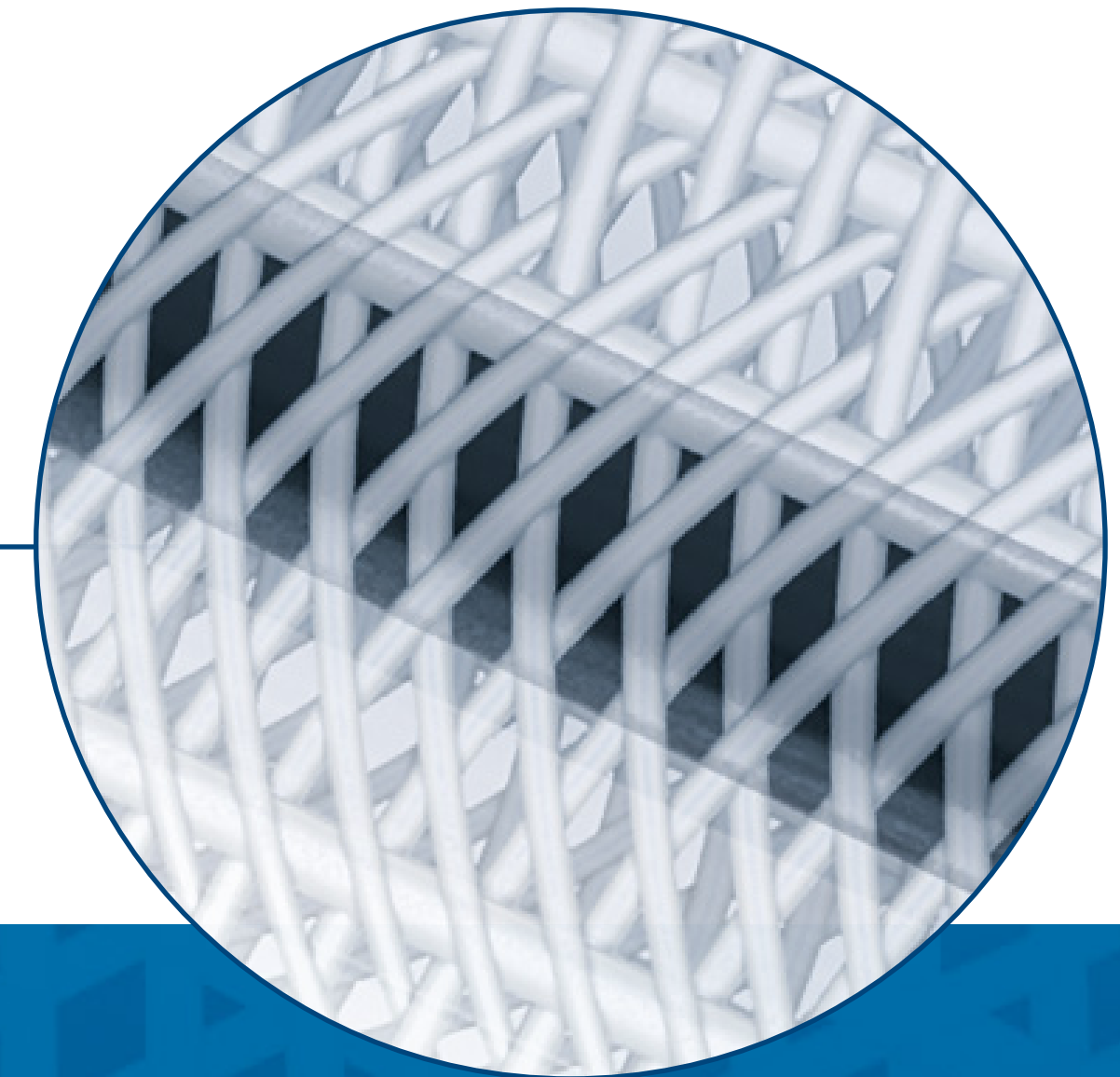


OPTIMIZED DELIVERABILITY

With an ideal blend of track and a low crossing profile, Athletis is designed to reach highly calcified, long, stenotic, or resistant fibrotic lesions, including those commonly found in AV access procedures.



ULTRA-LOW PROFILE

Athletis features a class-leading low lesion entry profile with an 0.042" flexible, tapered tip over an 0.035" guidewire.

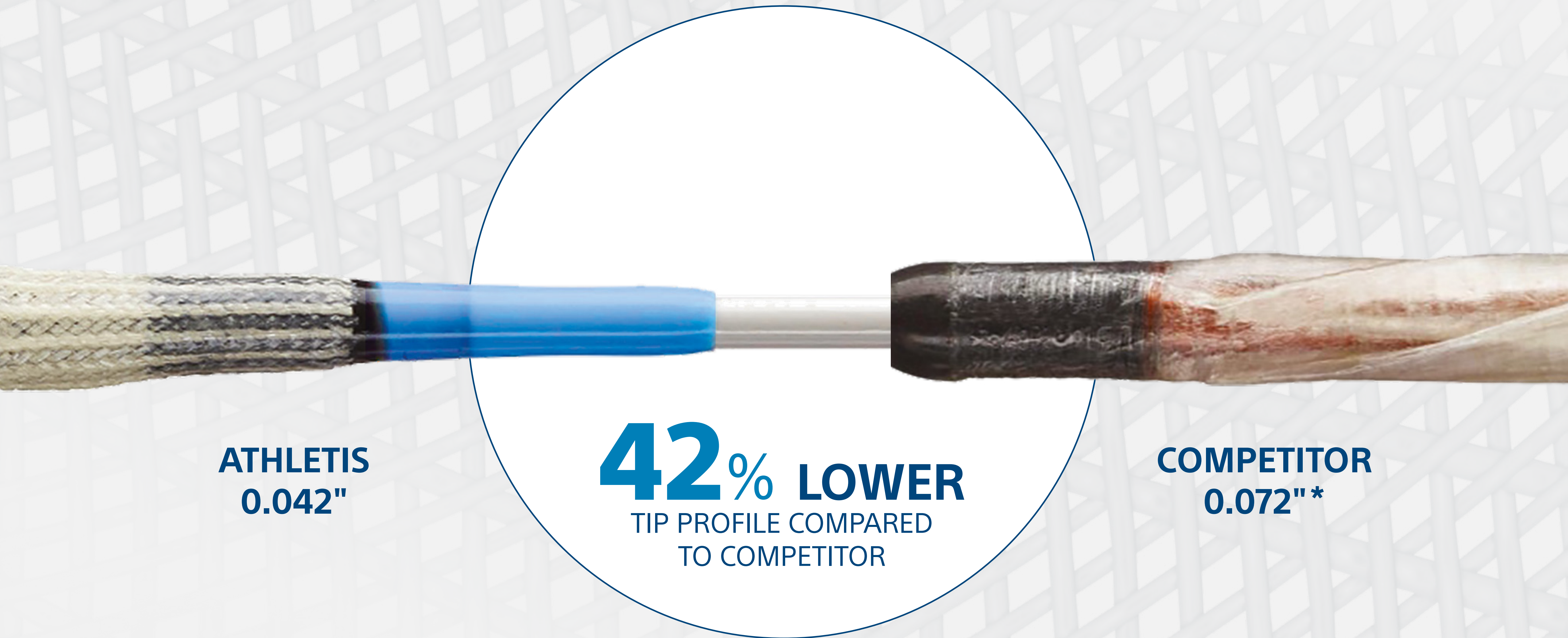

ATHLETIS
0.042"

42% LOWER
TIP PROFILE COMPARED
TO COMPETITOR

COMPETITOR
0.072"*



ATHLETIS ORDERING INFORMATION

UPN				Balloon Length (mm)					
				20	30	40	60	80	100
6F	7F	8F							
Balloon Diameter (mm)	4	Catheter Working Length (cm)	50	H74939347040250		H74939347040450	H74939347040650	H74939347040850	H74939347041050
			75	H74939347040270		H74939347040470	H74939347040670	H74939347040870	H74939347041070
	5		50	H74939347050250		H74939347050450	H74939347050650	H74939347050850	H74939347051050
			75	H74939347050270		H74939347050470	H74939347050670	H74939347050870	H74939347051070
			135			H74939347050410			
	6		50	H74939347060250		H74939347060450	H74939347060650	H74939347060850	H74939347061050
			75	H74939347060270		H74939347060470	H74939347060670	H74939347060870	H74939347061070
			135	H74939347060210		H74939347060410			
	7		50	H74939347070250		H74939347070450	H74939347070650	H74939347070850	H74939347071050
			75	H74939347070270		H74939347070470	H74939347070670	H74939347070870	H74939347071070
			135	H74939347070210		H74939347070410			
	8		50	H74939347080250	H74939347080350	H74939347080450	H74939347080650	H74939347080850	H74939347081050
			75	H74939347080270	H74939347080370	H74939347080470	H74939347080670	H74939347080870	H74939347081070
			135	H74939347080210		H74939347080410			
	9		50	H74939347090250		H74939347090450		H74939347090850	
			75	H74939347090270		H74939347090470		H74939347090870	
			135			H74939347090410			
	10		75	H74939347100270		H74939347100470		H74939347100870	
			135			H74939347100410			
	12		75	H74939347120270		H74939347120470	H74939347120670		
			135	H74939347120210		H74939347120410	H74939347120610		

INFLATION DEVICE ORDERING INFORMATION

Encore 40 is the recommended inflation device paired with Athletis. To order, use **UPN M001394472540**.

*Average measurements taken by Boston Scientific . Data on file (TR 92200478).
Bench test results may not necessarily be indicative of clinical performance.

All photographs taken by Boston Scientific

All cited trademarks are the property of their respective owners. CAUTION: The law restricts these devices to sale by or on the order of a physician. Indications, contraindications, warnings and instructions for use can be found in the product labeling supplied with each device. Information for the use only in countries with applicable health authority product registrations.

PI-856703-AA © 2020 Boston Scientific Corporation or its affiliates. All rights reserved.

